

Data Format Standardization and DICOM Integration for Hyperpolarized ^{13}C MRI

Ernesto Diaz¹, Renuka Sriram¹, Jeremy W. Gordon¹, Avantika Sinha¹, Xiaoxi Liu¹, Sule Sahin¹, Jason Crane¹, Marram P Olson¹, Hsin-Yu Chen¹, Jenna Bernard¹, Daniel B. Vigneron^{1,2}, Zhen Jane Wang¹, Duan Xu^{1,2}, Peder E. Z. Larson^{1,2}

Affiliations:

¹Department of Radiology and Biomedical Imaging, University of California – San Francisco, San Francisco, California, USA

²UC Berkeley-UCSF Graduate Program in Bioengineering, University of California, Berkeley and University of California, San Francisco, California, USA

Corresponding author address:

peder.larson@ucsf.edu, 1700 4th Street, Byers Hall Room 102C, San Francisco, CA 94143

Submitted to the *Journal of Imaging Informatics in Medicine* as a *Technical Note*

Keywords: Hyperpolarized ^{13}C MRI, metabolic imaging, DICOM format, Experiment metadata

Statements and Declarations

Funding

This work was supported by NIH grants P41EB013598, U24CA253377, R01CA262630, and R01CA249909. Authors PEZL, DBV, and JWG received research support from GE Healthcare.

Competing Interests

Financial interests: Authors ED, RS, AS, XL, JB, ZJW, and XU declare they have no financial interests. Authors PEZ, DBV, and JWG received research funding from GE Healthcare

Author Contributions

All authors contributed to the study conception and design. DICOM tools were created by Ernesto Diaz, Jeremy Gordon, Jason Crane, Marram P Olson, and Xiaoxi Liu. Data collection and analysis were performed by Ernesto Diaz. Creation of the proposed parameters was performed by Peder Larson. The first draft of the manuscript was written by Ernesto Diaz and Peder Larson, and all authors commented on previous versions of the manuscript. All authors read and approved the final manuscript.

Ethics approval

This manuscript presents technical methods but does not include any in vivo study results.

Consent to participate

Not applicable

Consent to publish

Not applicable

Abstract

Hyperpolarized (HP) ^{13}C MRI has shown promise as a valuable modality for *in vivo* measurements of metabolism and is currently in human trials at 15 research sites worldwide. With this growth it is important to adopt standardized data storage practices as it will allow sites to meaningfully compare data.

In this paper we (1) describe data that we believe should be stored and (2) demonstrate pipelines and methods that utilize the Digital Imaging and Communications in Medicine (DICOM) standard. This includes proposing a set of minimum set of information that is specific to HP ^{13}C MRI studies. We then show where the majority of these can be fit into existing DICOM Attributes, primarily via the “Contrast/Bolus” module.

We also demonstrate pipelines for utilizing DICOM for HP ^{13}C MRI. DICOM is the most common standard for clinical medical image storage and provides the flexibility to accommodate the unique aspects of HP ^{13}C MRI, including the HP agent information but also spectroscopic and metabolite dimensions. The pipelines shown include creating DICOM objects for studies on human and animal imaging systems with various pulse sequences. We also show a python-based method to efficiently modify DICOM objects to incorporate the unique HP ^{13}C MRI information that is not captured by existing pipelines. Moreover, we propose best practices for HP ^{13}C MRI data storage that will support future multi-site trials, research studies and technical developments of this imaging technique.

Introduction

Hyperpolarized (HP) ^{13}C MRI has shown promise as a valuable modality particularly for *in vivo* measurements of metabolism¹, and is currently in human trials at 15 research sites worldwide with over 75 human research publications². It is based on intravenous injection of a hyperpolarized contrast agent that has been enriched with ^{13}C . The most widely used agent is ^{13}C -pyruvate for measuring metabolism, but there are also a range of other promising agents undergoing clinical translation including ^{13}C -urea (perfusion)³, ^{13}C -alpha-ketoglutarate (metabolism)⁴, and ^{13}C -fumarate (necrosis)⁵. To date, human studies have been performed to characterize normal metabolism, and clinical studies have been performed to study in altered metabolism in prostate cancer, brain tumors, kidney tumors, pancreatic cancer, metastatic disease, liver disease, and heart disease².

As HP ^{13}C MRI continues to grow, it is important to adopt standardized data storage practices, allowing sites to compare data. The Digital Imaging and Communications in Medicine (DICOM) format is an attractive standard format as it is commonly used in digital imaging, medical imaging, and communications. The DICOM file type consists of the image and the metadata of the image packed into a single file. The information in the metadata is organized as a constant and is standardized by a series of data elements. By extracting the data elements, we can access important information regarding the patient demographics and study parameters that are crucial for study interpretation. The DICOM standard is actively supported and updated. It also supports the interchange of information between computer systems such as picture archiving and communication system (PACS).

The goal of this work is to describe an approach for HP ^{13}C MRI data standardization utilizing the DICOM format. We first describe the data requirements, focusing on the unique aspects of HP ^{13}C MRI. We then describe current processing pipelines that support the creation of DICOM objects from multiple vendors and pulse sequences. Finally, we demonstrate how this information can be incorporated into DICOM format.

Hyperpolarized ^{13}C MRI Data

The data requirements for HP ^{13}C MRI are outlined briefly below. This includes performing imaging to spatially resolve signals, encoding of multiple metabolites to measure metabolic conversion, performing dynamic measurements to capture the rapid kinetics, and characteristics of the HP agent, all of which that will affect the results and analysis.

Imaging

HP ^{13}C experiments can be spatially resolved using MRI techniques, providing important localization to visualize metabolic activity within the body. Therefore, a format to store multi-dimensional data needs to be supported.

Metabolite encoding

HP ^{13}C MRI has the unique ability to encode signals from multiple metabolites, distinguishing the injected agent, or substrate, from resulting metabolic products that are converted *in vivo*. The most common combination of metabolites is $[1-^{13}\text{C}]$ pyruvate (substrate), $[1-^{13}\text{C}]$ lactate (product), $[1-^{13}\text{C}]$ alanine (product), and ^{13}C -bicarbonate (product). Thus, the data should support a metabolite encoding dimension. Metabolite encoding in the acquisition is typically done with MR spectroscopy methods, chemical shift encoding (e.g. IDEAL) methods, or metabolite-

specific imaging⁶. The data must support metabolite encoding through identification of the individual compounds and/or the resonant frequencies, or with a spectroscopy dimension.

Dynamic Measurements

In a HP ¹³C MRI experiment, rapid metabolite kinetics following injection, including perfusion and metabolic conversion, that are typically captured with time-resolved imaging. Thus, the data should support dynamic measurements and record timings of the measurements.

HP ¹³C Agents

To analyze HP ¹³C data, it is important to know the characteristics of the HP ¹³C agent, how it was prepared, and how it was used. Table 1 lists such relevant HP ¹³C agent information. Thus, the data should support capturing of this information.

HP ¹³ C Agent Information	Purpose/Description
HP ¹³C agent(s) injected	Describe the composition of the injected agents. Note that this can include simultaneous injection of multiple agents
Administration route of agent	How was the HP solution delivered
Total Injected Volume	How much HP solution was injected in total
Concentration of each HP agent	Concentration of each HP agent in HP solution
Start Time of injection	Describe the timing of the start of injection, important to know relative to the data acquisition timing
Total injection duration	Over how long was the HP solution injected
Polarization	What was the measured polarization of the HP agent(s). Typically back-calculated to the time of dissolution.
Polarization Measurement Timing	timing of the polarization measurement relative to the data acquisition
Agent Relaxation Rates (e.g. T1)	Expected or measured relaxation rates used to calculate polarization
Dissolution Timing	Timing of the dissolution of the HP agent, most important to know relative to the data acquisition

Table 1:Hyperpolarized ¹³C agent information that we typically keep track of during our studies.

DICOM Processing Pipelines

Our current processing pipelines support DICOM creation from multiple vendors and pulse sequences, listed in Table 2. The spectroscopy sequences use an open-source standards-based software framework(SIVIC) for DICOM creation which includes support for DICOM Spectroscopy (7,8).

Scanner Vendor	Sequences	Raw Data	Tools	Output
Bruker	Slab-selective MRS, chemical shift imaging (CSI) and Echo-planar spectroscopic imaging	MRS Raw File or Dad File	SIVIC MATLAB	MR Spectroscopy DICOMs
Bruker	Metabolite-Specific Imaging	FID	Bruker Paravision MATLAB	Metabolite image, AUC, kinetic rate DICOMs
GE Healthcare	Slab-selective MRS, chemical shift imaging (CSI) and Echo-planar spectroscopic imaging	Pfile or ScanArchive	SIVIC	MR Spectroscopy DICOMs
GE Healthcare	Metabolite-specific EPI	Pfile or ScanArchive	GE Orchestra Toolbox MATLAB	Metabolite images, AUC ratio and kinetic rate DICOMs
RTHawk Research (Vista.ai)	Metabolite-specific spiral and bSSFP	RTHawk raw data Real-time reconstructed DICOM	MATLAB to improve data	Metabolite image DICOMs

Table 2: Currently supported sequences with DICOM creation pipelines in use at our institution in alphabetical order.

MR Spectroscopy and Spectroscopic Imaging with SIVIC

MR spectroscopy and spectroscopic imaging (MRS/I) for HP ^{13}C MRI includes data sampling at various time delays to resolve a spectrum and separate ^{13}C -labeled compounds. For processing of MRS/I data, we utilize the libraries in SIVIC for data processing and DICOM creation, including support of the DICOM MRS standard⁹. This software can reads in multiple vendor MRS raw data files (Bruker, GE, Philips, Siemens, Varian), and it for HP ^{13}C MRI studies on GE Healthcare and Bruker MRI scanners.

1. Reads in raw data
2. Performs MR spectroscopy or MR spectroscopic imaging reconstructions

3. Writes to DICOM MRS file
4. Optional – extract peak amplitudes of areas from spectra and create individual metabolite DICOM files
5. Optional - write raw data or partially processed data to DICOM
6. Optional - secondary capture of MRS and anatomical visualization to DICOM

GE Metabolite-specific Imaging

Metabolite-specific imaging for HP ^{13}C MRI uses a spectral-spatial RF pulse to selectively excite a single metabolite¹⁰ and is followed by fast imaging readout such as echo-planar imaging (EPI) or spirals. We have implemented EPI-based metabolite-specific imaging for GE Healthcare MRI scanners along with a reconstruction pipeline based on the GE Orchestra reconstruction toolbox¹¹. This leverages the EPI processing and DICOM creation functionality that is utilized by GE product pulse sequences.

1. Reads in raw data (Pfile or ScanArchive)
2. Perform EPI reconstruction including Nyquist ghost correction
3. Perform HP-specific coil combination¹²
4. Create area-under-curve (AUC) and kinetic rate maps
5. Write metabolite images and parameter maps to DICOM

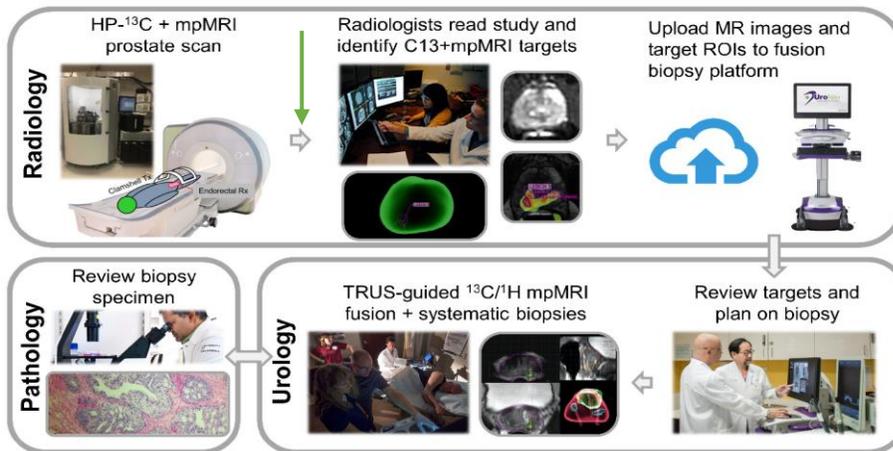


Figure 1: Example of how HP ^{13}C MRI can be used in a workflow for improved targeting of biopsies in prostate cancer¹¹. Using DICOM has been critical in the success of this workflow, as this allows for integration into PACS as well as the fusion biopsy platform, with radiologists and urologists easily able to review the images and resulting annotations outlining the biopsy targets. Figure reproduced with permission from¹¹.

RTHawk Metabolite-specific Imaging

We have implemented metabolite-specific imaging with EPI and spirals, as well as metabolite-specific balanced steady-state free-precession (bSSFP) sequences¹³⁻¹⁵ on the RTHawk Research platform. This is a vendor-neutral platform that allows for advanced control of the

MRI system including real-time reconstructions and feedback as well as pulse sequence programming.

1. We utilize real-time reconstruction built into RTHawk during the scanning process, resulting in individual metabolite DICOM images at each timepoint. (Note: This uses the sum of squares coil combination technique.)
2. The DICOMs Images are exported when the data is transferred from the RTHawk workstation to the local network.
3. A second, improved image reconstruction is performed using raw complex-valued data and a more advanced coil combination method to enhance the data quality⁹.
4. To integrate the improved data into DICOM, a MATLAB script is employed that reads the metadata from the real-time reconstructed DICOMs Images, replaces the original voxel data with the results of the second, optimized reconstruction, and creates new DICOMs Images with appropriate new unique identifiers (UIDs).

Bruker Metabolite-specific Imaging

Similar to the GE metabolite-specific imaging sequence, it uses a spectral-spatial RF pulse to excite a single metabolite¹⁰, and is followed by a fast echo-planar imaging (EPI) readout¹⁶.

1. Bruker scanner has its own built-in EPI reconstruction and an option to write reconstructed metabolite maps to DICOM files.
2. A MATLAB script is used to read in the DICOM files. From that data, we generate AUC maps and kinetic rate maps, and write back to DICOM files.
3. Optional: Reads in Bruker raw data using a MATLAB script, use a custom reconstruction script and write to DICOM.

DICOM HP Agent Integration

Many of the HP agent metadata that are desirable for analyses (Table 1) can be integrated into existing DICOM attributes, particularly within the “Contrast/Bolus” module. However, the above DICOM pipelines by default do not include support for these metadata. To address this, we built a tool to modify HP ¹³C DICOM objects based on user input to include this HP agent metadata.

Currently, the HP ¹³C images created by MRI scanners or in-house reconstruction pipelines do not utilize the Contrast/Bolus Attributes. We seek to leverage the existing Contrast/Bolus Attributes in DICOM to store information to identify the study and experimental parameters that are valuable for HP ¹³C data analysis by researchers and clinicians. Table 3 lists the attributes we are adding to our images.

TAG	VR	Attribute Name	¹ H MRI Example Values	HP ¹³ C Description	HP ¹³ C Example Values
(0018,0010)	LO	Contrast/Bolus Agent	CC MAGNEVIST	HP ¹³ C agent(s) injected	HYPERPOLARIZED [1- 13C]PYRUVATE
			ORAL & OMNIPAQUE		HYPERPOLARIZED [2- 13C]PYRUVATE

					HYPERPOLARIZED [1- ¹³ C]PYRUVATE + [13C,15N]-UREA
(0018,1040)	LO	Contrast/Bolus Route		Administration route of contrast agent	IV
(0018,1041)	DS	Contrast/Bolus Volume		Injected volume (mL)	
(0018,1042)	TM	Contrast/Bolus Start Time		Define start of injection relative to data acquisition (Not DICOM compliant) Start Time in HHMMSS.FFFFFFFF (DICOM compliant)	
(0018,1047)	DS	Contrast/Bolus Flow Duration		Total injection duration [s]	
(0018,1048)	CS	Contrast/Bolus Ingredient	IODINE BARIUM GADOLINIUM CARBON DIOXIDE	¹³ C enriched compounds	[1- ¹³ C]PYRUVATE [2- ¹³ C]PYRUVATE [¹³ C, ¹⁵ N]UREA [1- ¹³ C]ALPHA-KETOGLUTARATE
(0018,1049)	DS	Contrast/Bolus Ingredient Concentration		Concentration of each ¹³ C compound [mg/mL]	
(0400,0550)	SQ	Modified Attributes Sequence		Store prior versions of an attributes that were removed or modified	
(0400,0562)	DT	Attribute Modification DateTime		Tracked Date and Time when attributes were modified	
(0400,0563)	LO	Modifying System		Describes what modified the attributes	SIVIC-HP_agent_DICOM_tool.py

Table 3: DICOM attributes proposed for addition to HP ¹³C MRI study data. This captures the majority but not all the desired parameters listed in Table 1. Note that the “Contrast/Bolus Ingredient” attributes listed can include more than one entry to support multiple simultaneously polarized HP ¹³C agents¹⁷ such as pyruvate/urea^{3,14}. The proposed ¹³C naming conventions are in a style that is consistent with the style used for PET in DICOM.

In DICOM, each attribute or data element is identified by a unique tag (TAG), which is comprised of a group number and an element number. The value representation (VR) is a code that indicates the data type used to encode the value(s) of the attribute. Additionally, the value multiplicity (VM) indicates how many values can be present in the attribute. These three

components, work together to provide a comprehensive understanding of the attribute and its role within the DICOM file, allowing for accurate and efficient manipulation and interpretation of the medical imaging data.

We reviewed the DICOM standards definition and found the listed attributes in Table 3 can capture much of the unique HP agent metadata. This includes ingredients, injection timing, and dosing characteristics. These can all be added in a DICOM compliant fashion.

One important feature captured in this proposal is the inclusion of multiple ingredients, which can capture dual-agent (“co-polarized”) or multi-agent HP studies^{3,14,17}. In these studies, multiple ¹³C-enriched compounds are simultaneously hyperpolarized and simultaneously injected. The “Contrast/Bolus Agent Sequence”, which can capture use of multiple agents that affect the data in a DICOM object, was not used since when multiple HP compounds are desired, they are combined into a single agent with multiple ingredients. The use of multiple agents is not performed to our knowledge, as it would require rapid dissolution of multiple separate samples following hyperpolarized to inject both before the signal decays.

We also proposed standardized ingredient names that are consistent with the style used in PET. In these, the enriched isotope(s) are denoted with “^” around the nucleus, e.g. “¹³C”. We also use a standard chemical naming convention to identify the site of enrichment within the molecule.

We used the PyDicom python toolbox to access the data in DICOM files (Figure 2). PyDicom is a general-purpose DICOM framework whose purpose is to reads and write DICOM attributes. It can also add, delete, and modify any attributes to DICOM objects.

```
reading  ───> path = '/Users/ernestodiaz/Desktop/DICOM/AUC_lactate/Image_001.dcm'
adding   ───> ds = pydicom.dcmread(path)
writing  ───> ds.add_new([0x0010,0x0010], 'PN', "Ernesto")
           ds.save_as(path)
           print(ds)
           (0010, 0010) Patient's Name                PN: 'Ernesto' ← results
```

Figure 2: This is an example of the PyDicom functions used. It reads a file, then adds an attribute and saves the file.

The basic workflow of the python script created using the packages PyDicom and OS is as follows: User would run the script in the terminal with a parameter of a directory path for DICOMS that need to be modified. The prompts would pop up one by one and the user can write its own inputs or press enter for default values as shown in Figure 3. Once the user has entered the prompts, it will loop through the folder and add all the attributes to each DICOM file.

We have also aimed to capture the modifications made to the DICOM file when the tool is used. For this, we used the Modified Attributes Sequence to describe what attributes we modified. For Modifying System, it shows the name of this tool we used to change the attributes. We used DateTime within the Modified Attributes Sequence to keep track of the date when we modified

the attributes, and this shows other researchers when we modified them. This maybe a non-standard use of these attributes.

This tool has been demonstrated on DICOM objects from a metabolite-specific EPI sequence¹⁸ on a GE 3T scanner, reconstructed with a custom implementation based on the GE Orchestra package, DICOM from Bruker 2D CSI scan generated by SIVIC, and a DICOM from a GE echo-planar spectroscopic imaging (EPSI) sequence generated by SIVIC.

The script is available within SIVIC at https://github.com/SIVICLab/sivic/tree/master/applications/dicom_tools.

Figure 4: The prompts questions that gives option for default values or user input. As well, showing the attributes that has been modified.

```
----Welcome----
Enter A Contrast Agent (Press Enter For Default: HYPERPOLARIZED [1-13C]PYRUVATE):
The Current Administration Route of Contrast Agent is Oral & IV. Would you like to overwrite it? Type Y or N:
Enter The Volume Of The Contrast Agent Injected In mL (Press Enter For Default: 0.35):
Enter The Total Injection Duration In Seconds (Press Enter For Default: 12):
Enter The Delay Time From Start Of Injection To Start Of Acquisition In Seconds (Press Enter For Default: 0):
Enter A Contrast Ingredient (Press Enter For Default: [1-^13^C]PYRUVATE):
Enter The Molarity Of The Contrast Ingredient in mM (Press Enter For Default: 80):
Enter The Molar Mass Of The Contrast Ingredient in g/mol (Press Enter For Default: 88.06):
Would You Like To Add Another Ingredient? Type Y or N:
IM-0001-0241.dcm
These are the attributes that you modified
[(0400, 0562) Attribute Modification DateTime      DT: '2024031623'
(0400, 0563) Modifying System                    L0: 'SIVIC-HP_Agent_DICOM_Tool_Test.py']
```

HP ¹³C MRI Data Storage Best Practices

We have also created a proposed set of Data Storage Best Practices. These were based upon a survey of research groups at our institution, performing both human and preclinical HP ¹³C imaging studies at our institution. The proposed set of Data Storage Best Practices aim to ensure that the stored data completely describes the study as is necessary for subsequent processing. This should also be done so that future researchers can easily retrieve and analyze the data without the need to search multiple sources of information.

Storage of essential study data should be done in a centralized location with high-quality back up, using a single folder per study, including:

- HP spectroscopy/image data (DICOM recommended)
- HP raw data
- ¹H image data (DICOM recommended)
- Copy of any processing scripts
- Metadata – any other information about the study. This can either be directly noted in the data folder, or in a separate database/spreadsheet. At a minimum this should include
 - Image acquisition data (ideally in DICOM files):
 - Flip angle information, especially for variable flip angle or metabolite-specific flip schedules
 - Injection timing characteristics – length of injection, timing of imaging relative to injection(See Table 1)
 - Dose characteristics – volume injected, concentration of agent, subject weight, polarization.(See Table 1)

- Other Study data - directly recorded in the study folder, or in a separate database/spreadsheet.
- Study “Key” needed in all data locations to link to other metadata sources.
 - Option 1: Unique ID for each study
 - Option 2: Unique ID for each subject + Unique ID for each study

Discussion and Unmet Needs

Using the DICOM standard has significant advantages for HP ^{13}C MRI data. Given this technology is being translated into clinical studies, DICOM will allow for seamless integration into existing image viewing and analysis platforms, ultimately supporting use in clinical workflows. The use of standardized DICOM objects would also be valuable in the multi-site clinical trials setting, which are likely to begin within the next several years.

We also believe that the DICOM format will be valuable for the preclinical studies given the flexibility of the format. A pre-clinical DICOM standard has recently been proposed¹⁹, and a recent assessment preclinical imaging metadata need has also advocated for integration into DICOM²⁰. While adoption is less consistent in preclinical imaging systems, we have demonstrated a pipeline for our Bruker preclinical MRI scanner.

While the Contrast/Bolus module in DICOM supported most HP agent parameters we wanted to save, there were several parameters without an existing attribute. These include the polarization measurements, T1 relaxation rate, and dissolution timing. These could be stored in private data elements. This could also justify the need for a HP Agent MRI Module in DICOM. This could largely follow the examples from PET, which also must capture similar information: polarization is analogous to specific activity, relaxation rates are analogous to half-life, and dissolution timing is analogous to dose creation timing. The PET Isotope Module (C.8.9.2) could serve as a good model.

Several other current features of the DICOM standard may also be relevant for HP MRI. The Enhanced MR Information Object Definitions (IODs) are attractive because they allow for more MRI specific data sources such as raw k-space and spectroscopy data and can be used for “raw” images (e.g. dynamic metabolite images) as well as derived parameter maps (e.g. kinetic rate maps). The Enhanced Contrast/Bolus Module may also be relevant to store additional HP agent parameters.

Hyperpolarization is also used in HP ^{129}Xe MRI, a modality that unique images pulmonary ventilation, gas exchange, and terminal airway morphology²² and received FDA approval in 2021. There are many shared aspects of both HP MRI methods, including the hyperpolarization of an agent, administration of this agent (inhaled for HP ^{129}Xe), relatively rapid imaging of the agent, and measuring agent kinetics. With this in mind, we believe it would be beneficial to coordinate efforts across the HP modalities when developing DICOM integration methods and specifications of metadata.

HP ^{13}C MRS and MRSI acquisition methods are more challenging to use in a standardized fashion compared to imaging-based methods (e.g. metabolite-specific imaging). This is because the adoption of DICOM MRS by vendors is rather limited, leading to the community development of tools such as SIVIC⁸. If MRS and MRSI methods continue to show value, it

would be important for researchers to work closely with vendor partners and ideally have the DICOM MRS standard integrated into the scanners.

Ultimately, we believe it is important that both a standardized and essential set of DICOM attributes are recorded for HP ¹³C MRI studies. This should be standardized so as not to introduce confusion, for example due to different naming conventions or use of different DICOM attributes. This should also include an essential set of attributes to allow for robust and reproducible analyses of HP ¹³C MRI data. The final definitions are beyond the scope of this paper, but we think would be best addressed by a consensus of HP ¹³C MRI researchers in future discussions as the community recently formed the HP 13C Consensus Group². Following consensus, any changes to the DICOM standard should then be proposed by drafting a Correction Proposal (CP) and work with the DICOM MR Working Group (WG 16) to submit it to the base standards WG (WG 6).

Conclusion

We have proposed a set of attributes to store HP ¹³C MRI metadata, particularly regarding the HP agent, and implemented via PyDicom a program to add these a set of these attributes in a DICOM-compliant fashion based on user input. This was shown to work in line with a variety of DICOM creation pipelines that are used for HP ¹³C data. In addition, we have proposed a broader set of data storage best practices. Providing further data standardization will advance this technology by allowing sites to seamlessly share and compare data, including in the context of multi-site trials, and by using DICOM HP ¹³C MRI can be easily integrated into other clinical and research workflows.

References

1. Wang ZJ, Ohliger MA, Larson PEZ, et al. Hyperpolarized ^{13}C MRI: State of the Art and Future Directions. *Radiology*. Published online March 5, 2019:182391. doi:10.1148/radiol.2019182391
2. Larson PE, Bernard JM, Bankson JA, et al. Current Methods for Hyperpolarized [1- ^{13}C]pyruvate MRI Human Studies. Published online September 7, 2023. doi:10.48550/arXiv.2309.04040
3. Qin H, Tang S, Riselli AM, et al. Clinical translation of hyperpolarized ^{13}C pyruvate and urea MRI for simultaneous metabolic and perfusion imaging. *Magnetic Resonance in Medicine*. n/a(n/a). doi:10.1002/mrm.28965
4. Chaumeil MM, Larson PEZ, Yoshihara HAI, et al. Non-invasive in vivo assessment of IDH1 mutational status in glioma. *Nat Commun*. 2013;4:2429. doi:10.1038/ncomms3429
5. Clatworthy MR, Kettunen MI, Hu DE, et al. Magnetic resonance imaging with hyperpolarized [1,4-(^{13}C) $_2$]fumarate allows detection of early renal acute tubular necrosis. *Proc Natl Acad Sci U S A*. 2012;109(33):13374-13379. doi:10.1073/pnas.1205539109
6. Gordon JW, Chen H, Dwork N, Tang S, Larson PEZ. Fast Imaging for Hyperpolarized MR Metabolic Imaging. *Magnetic Resonance Imaging*. 2021;53(3):686-702. doi:10.1002/jmri.27070
7. Crane JC, Olson MP, Nelson SJ. SIVIC: Open-Source, Standards-Based Software for DICOM MR Spectroscopy Workflows. *International Journal of Biomedical Imaging*. 2013;2013:e169526. doi:10.1155/2013/169526
8. SIVIC. doi:10.5281/zenodo.4777197
9. Crane JC, Gordon JW, Chen HY, et al. Hyperpolarized ^{13}C MRI data acquisition and analysis in prostate and brain at University of California, San Francisco. *NMR in Biomedicine*. 2020;n/a(n/a):e4280. doi:10.1002/nbm.4280
10. Cunningham CH, Chen AP, Lustig M, et al. Pulse sequence for dynamic volumetric imaging of hyperpolarized metabolic products. *J Magn Reson*. 2008;193(1):139-146.
11. Chen HY, Bok RA, Cooperberg MR, et al. Improving multiparametric MR-transrectal ultrasound guided fusion prostate biopsies with hyperpolarized ^{13}C pyruvate metabolic imaging: A technical development study. *Magn Reson Med*. 2022;88(6):2609-2620. doi:10.1002/mrm.29399
12. Zhu Z, Zhu X, Ohliger MA, et al. Coil combination methods for multi-channel hyperpolarized ^{13}C imaging data from human studies. *Journal of Magnetic Resonance*. 2019;301:73-79. doi:10.1016/j.jmr.2019.01.015
13. Tang S, Bok R, Qin H, et al. A metabolite-specific 3D stack-of-spiral bSSFP sequence for improved lactate imaging in hyperpolarized [1- ^{13}C]pyruvate studies on a 3T clinical scanner. *Magnetic Resonance in Medicine*. n/a(n/a). doi:10.1002/mrm.28204
14. Liu X, Tang S, Mu C, et al. Development of specialized magnetic resonance acquisition techniques for human hyperpolarized [^{13}C , $^{15}\text{N}_2$]urea + [1- ^{13}C]pyruvate simultaneous

perfusion and metabolic imaging. *Magn Reson Med.* 2022;88(3):1039-1054.
doi:10.1002/mrm.29266

15. Liu X, Tang S, Cui D, et al. A metabolite specific 3D stack-of-spirals bSSFP sequence for improved bicarbonate imaging in hyperpolarized [1-¹³C]Pyruvate MRI. *J Magn Reson.* 2023;353:107518. doi:10.1016/j.jmr.2023.107518

16. Sahin SI, Ji X, Agarwal S, et al. Metabolite-Specific Echo Planar Imaging for Preclinical Studies with Hyperpolarized ¹³C-Pyruvate MRI. *Tomography.* 2023;9(2):736-749.
doi:10.3390/tomography9020059

17. Wilson DM, Keshari KR, Larson PEZ, et al. Multi-compound Polarization by DNP Allows Simultaneous Assessment of Multiple Enzymatic Activities In Vivo. *J Magn Reson.* 2010;205(1):141-147.

18. Gordon JW, Vigneron DB, Larson PEZ. Development of a symmetric echo planar imaging framework for clinical translation of rapid dynamic hyperpolarized (¹³C) imaging. *Magn Reson Med.* Published online February 2016. doi:10.1002/mrm.26123

19. Kalen JD, Clunie DA, Liu Y, Tatum JL, Jacobs PM, Kirby J, Freymann JB, Wagner U, Smith KE, Suloway C, et al. Design and Implementation of the Pre-Clinical DICOM Standard in Multi-Cohort Murine Studies. *Tomography.* 2021; 7(1):1-9.
<https://doi.org/10.3390/tomography7010001>

20. Moore SM, Quirk JD, Lassiter AW, Laforest R, Ayers GD, Badea CT, Fedorov AY, Kinahan PE, Holbrook M, Larson PEZ, et al. Co-Clinical Imaging Metadata Information (CIMI) for Cancer Research to Promote Open Science, Standardization, and Reproducibility in Preclinical Imaging. *Tomography.* 2023; 9(3):995-1009. <https://doi.org/10.3390/tomography9030081>

21. Niedbalski PJ, Hall CS, Castro M, et al. Protocols for multi-site trials using hyperpolarized ¹²⁹Xe MRI for imaging of ventilation, alveolar-airspace size, and gas exchange: A position paper from the ¹²⁹Xe MRI clinical trials consortium. *Magn Reson Med.* 2021; 86: 2966–2986.
<https://doi.org/10.1002/mrm.28985>